










<p><u>Big Button Telephone</u></p> <p><u>Standalone item</u> (does not talk to other TEC items)</p> <p>CES71009-0052</p> 	<p>Large, buttoned telephone. Standalone device.</p> <p>Makes it easier to see & manipulate to make calls.</p> <p>Sound amplifier makes telephone louder.</p> <p>Has answerphone function.</p> <p>Where you can find it on the Livity Life ordering system;</p> <p>Category: Dexterity Aid</p>
<p><u>Smoke alarm</u></p> <p>Linked to monitoring centre / pager.</p> <p>CES65736-0451</p> 	<p>Detects smouldering and open fires.</p> <p>Will notify Person of possible fire by sounding an audio alarm.</p> <p>Wireless smoke alarm linked Kit – links to lifeline unit (monitoring centre) or carer pager (range 500metres) which need to be ordered separately (Lifeline CES65742-0451 or Buzz CES65748-1645).</p> <p>Where you can find it on the Livity Life ordering system;</p> <p>Category: Environmental Sensors</p>
<p><u>Heat Detector</u></p> <p>Linked to monitoring centre / pager.</p> <p>CES65786-0451</p> 	<p>Detects extreme rapid rise temperature.</p> <p>Will notify person of possible fire by sounding an audio alarm.</p> <p>Wireless Heat Detector linked Kit – links to lifeline unit (monitoring centre) or carer pager (range 500metres) which need to be ordered (Lifeline CES CES65742-0451 / Buzz CES65748-1645).</p> <p>Where you can find it on the Livity Life ordering system;</p> <p>Category: Environmental Sensors</p>



<p><u>Carbon Monoxide detector</u></p> <p>Linked to monitoring centre / pager.</p> <p>CES65933-1416</p> 	<p>Detector will trigger if high levels of carbon monoxide are present.</p> <p>Wireless Carbon Monoxide alarm linked Kit – links to lifeline unit (monitoring centre) or carer pager (range 500metres) which need to be ordered (Lifeline CES CES65742-0451 / Buzz CES65748-1645).</p> <p>Where you can find it on the Livity Life ordering system;</p> <p>Category: Environmental Sensors</p>
<p><u>Flood & Water Detector</u></p> <p>Linked to monitoring centre / pager.</p> <p>CES65838-045</p> 	<p>Siren sounds when flood is detected.</p> <p>Can be prescribed as a standalone item – siren will alert other household members.</p> <p>Can be linked to lifeline unit (monitoring centre) or carer pager (range 500metres) which need to be ordered (Lifeline CES CES65742-0451 / Buzz CES65748-1645).</p> <p>Where you can find it on the Livity Life ordering system;</p> <p>Category: Environmental Sensors</p>
<p><u>Epilepsy Monitor -Companion Mini</u></p> <p>CES65788-0451</p> 	<p>Bed sensor with Pager – which alerts when tonic-clonic seizures are detected.</p> <p>It is placed beneath the mattress on either a divan or slatted bed.</p> <p>Works on 'movement'.</p> <p>Pager included with kit.</p> <p>Pager will alert Carer if sensor detects seizure activity</p> <p>Where you can find it on the Livity Life ordering system;</p> <p>Category: Epilepsy</p>

<p><u>Movement Sensitive cordless Light</u></p> <p><u>Standalone item</u> (does not talk to other TEC items)</p> <p>CES65966-0162</p> 	<p>Illuminates indoor and outdoor spaces.</p> <p>Light is activated when motion / movement is detected.</p> <p>Light is wireless & LED (bright)</p> <p>Motion / light activated can be detected up to 4 meters away.</p> <p>Light will remain on for as long as there is a presence in the beam.</p> <p>Automatically switches off 20 seconds after movement ceases.</p> <p>Where you can find it on the Livity Life ordering system;</p> <p>Category: Monitoring Movement in the home</p>
<p><u>Movement Detector - PIR Transmitter & Pager Kit-</u></p> <p><u>Standalone item</u> (does not talk to other TEC items)</p> <p>CES65632-0162</p> 	<p>PIR (beam) detects and alerts upon movement – detector talks to pager (both detector and pager come with this kit)</p> <p>Can be strategically placed depending on need; front door, bed, chair etc.</p> <p>Notifies Carer, via pager, when the person is moving/independently transferring etc when detector is triggered.</p> <p>Pager has a range of approx. 100metres.</p> <p>Movement detector can be wall mounted or comes with a stand.</p> <p>Where you can find it on the Livity Life ordering system;</p> <p>Category: Monitoring Movement in the home</p>

<p><u>Vibby Falls Detector</u></p> <p><u>CES65738-0451</u></p> 	<p>Automatically detects falls based upon movements – Wearer does not need to press a button to call for help, but this option is also available (they can push a button).</p> <p>Is worn around wrist.</p> <p>Cannot be worn around the neck.</p> <p>Water resistant.</p> <p>This is not a standalone item. Needs to be linked to either; lifeline unit (monitoring centre) or carer pager (range 500metres). which need to be ordered (Lifeline CES CES65742-0451 Buzz CES65748-1645).</p> <p>Where you can find it on the Livity Life ordering system;</p> <p>Category: Monitoring Movement in the home</p>
<p><u>Chair Leaving Alarm Kit with Voice Recorder</u></p> <p><u>CES65631-0162</u></p> 	<p>Kit comes with sensor mat for chair, pager and voice transmitter.</p> <p>Upon movement the voice transmitter will offer helpful individualised prompts – such as 'Dad remember to use your walking frame' whilst an alert is sent to pager to notify Carer.</p> <p>Pager range is approx. 100 metres.</p> <p>Item can be prescribed as standalone or used with Bed Leaving Alarm kit (CES65630-0162) and set up so two sensors will talk to one pager.</p> <p>Item does not talk to Monitoring Centre.</p> <p>Where you can find it on the Livity Life ordering system;</p> <p>Category: Monitoring Movement in the home</p>

<p><u>Bed Alarm Kit with Voice Recorder</u></p> <p><u>CES65630-0162</u></p>  <p>CT-3 BED SENSOR 75cm x 45cm</p> <p>BED ALARM KIT</p>	<p>Kit comes with sensor mat for bed, pager and voice transmitter.</p> <p>Upon movement the voice transmitter will offer helpful individualised prompts – such as ‘Dad don’t get out of bed, I’m on my way’ whilst an alert is sent to pager.</p> <p>Pager range is approx. 100 metres.</p> <p>Item can be prescribed as standalone or used with Chair Leaving Alarm kit (CES65631-0162).</p> <p>Item does not talk to Monitoring Centre</p> <p>Where you can find it on the Livity Life ordering system;</p> <p>Category: Monitoring Movement in the home</p>
<p><u>Bed Leaving Sensor – Linked to monitoring Centre.</u></p> <p>This is not a standalone item.</p> 	<p>Bed sensor works on pressure and will alert when pressure is lifted from the sensor – when the person transfers from the bed.</p> <p>Prescriber can ask Livity Life to set timings on this sensor – for instance, alert to sound when person has been out of bed for 20 minutes between 11pm & 7am.</p> <p>The alert can be set to trigger the pager (CES65748-1645) or if the pager is out of range / or monitoring centre .</p> <p>If monitoring centre support is required a lifeline (CES CES65742-0451) will be required.</p> <p>Where you can find it on the Livity Life ordering system;</p> <p>Category: Monitoring Movement in the home</p>

<p><u>Jelly Bean Neat 4066</u></p> <p>This is not a standalone item. CES74250-0531</p> 	<p>A button that triggers an alert / alarm. Button comes with coloured overlays (Red, Blue, Yellow, Green).</p> <p>Large Button – helpful for those with reduced hand function / poor dexterity.</p> <p>This is not a standalone item. Needs to be linked to either; lifeline unit (monitoring centre) or carer pager (range 500metres) which must be ordered (Lifeline CES CES65742-0451 Buzz CES65748-1645).</p> <p>For monitoring purposes responders will be required who live within 30 mins of the property.</p> <p>Where you can find it on the Livity Life ordering system;</p> <p>Category: Lifeline Accessories</p>
<p><u>Telecare - Smile Pendant White kit</u></p> <p>This is not a standalone item. CES65785-0451</p> 	<p>Additional button to accompany lifeline kit that comes with one button (CES CES65742-0451).</p> <p>The Smile pendant can be worn on wrist or around neck.</p> <p>Pendant is waterproof.</p> <p>This is not a standalone item. Needs to be linked to either; lifeline unit (monitoring centre) or carer pager (range 500metres) which need to be ordered (Lifeline CES65742-0451 Buzz CES65748-1645).</p> <p>Where you can find it on the Livity Life ordering system;</p> <p>Category: Lifeline Accessories</p>

<p><u>Medication Dispenser:</u></p> <p><u>Standalone</u></p> <p>CES65746-0451</p> 	<p>Tablets are locked in a secure compartment and automatically released at pre-set times. Up to 6 daily audio & visual alarms (flashing LED on the lid). Buzzer (30 minutes) sounds continuously until pills are dispensed. The alarm is cancelled by tipping the product to dispense the pills. 28 tablet compartments 3 alarm tones with excellent level of volume.</p> <p>Where you can find it on the Livity Life ordering system;</p> <p>Category: Medication Dispenser</p>
<p><u>Medication Dispenser:</u></p> <p><u>Monitored</u></p> <p>CES65746-0451</p> 	<p>Automatically dispenses pills according to pre-programmable times/days.</p> <p>28 Pill slots.</p> <p>Pre-Programmable set times – Reminds Customer to take medication by audio alarm & flashing light.</p> <p>This item can be linked to monitoring centre / Carer informed if pills are not dispensed at pre-programmed time.</p> <p>To link order either; lifeline unit (monitoring centre) or carer pager (range 500metres) which need to be ordered (Lifeline CES CES65742-0451 Buzz CES65748-1645).</p> <p>Where you can find it on the Livity Life ordering system;</p> <p>Category: Medication Dispenser</p>

Lifeline

Monitored
CES65742-0451



Lifeline unit with one pendant Alarm.

Can be worn around the neck or wrist.

Pendant alarm is waterproof; can be fully submerged for approx 20 mins.

Additional pendant alarms can be prescribed if needed (see 'Help Button', 'Smile Pendant CES65785-0451).

When pendant alarm is triggered, the Person will be able to talk to monitoring centre &/or monitoring centre will be aware of the alert if the person is not able to communicate.

Additionally, this Lifeline can be linked to a pager CES65748-1645) Pager will be alerted if alarm is triggered – but only if within a 500metre distance.

Where you can find it on the Livity Life ordering system;

Category: Monitoring Movement in the home

Canary Care Passive Monitoring System



Passive monitoring system. Kit comes with a hub unit & 6 sensors; including 2 which can be used as door sensors.

Sensors monitor movement, temperature and light.

The hub unit plugs into a socket, centrally in the home. Prescriber to state where the sensors need to be placed.

Does not need a phone line or internet connection.

Prescriber will have access to online portal to view data and activity charts.

Based upon Prescribers assessment access to the online portal can be granted to appropriate family members / Carers.

Can be set up alerts are triggered when 'rules' are broken – such as if the front door opens at 2am an appropriate family member may be informed.

Prescriber to set appropriate rules dependant on client need. This item is for assessment purposes and is typically used for 2 weeks. If a longer-term need if identified the prescriber should discuss directly the TEC Team.

Prescribers of Canary Care can contact IntegratingTEC@centralbedfordshire.gov.uk for support if required.

Where you can find it on the Livity Life ordering system;

Category: Monitoring Movement in the home

Pager- Care Buzzbox - Pager

CES65748-1645



This is **not a standalone** item.



Pager for Carer / family member to use.



Pager can be linked to multiple sensors.




Pager has a range of approx. 500metres.

Where you can find it on the Livity Life ordering system;

Category: Monitoring Movement in the home

<p><u>Epilepsy Monitor -Companion Mini</u></p> <p><u>Item Code: CES65784-0451</u></p> 	<p>Compact, robust and portable, the Companion Mini supports the detection of tonic-clonic seizures.</p> <p>It is placed beneath the mattress on either a Divan or Slatted base bed.</p> <p>This sensor can differentiate between seizure activity and normal sleep activities before sending an alert.</p> <p>Where you can find it on the Livity Life ordering system;</p> <p>Category: Epilepsy</p>
<p><u>Memrabel 2 / Dementia Friendly Clock</u></p> <p><u>CES65668-0162</u></p> 	<p>Dementia friendly clock – large display showing time of day; 11am Morning. Supports those with can be disorientated with time, day &/or date.</p> <p>Can set up to host 20 'Reminders'.</p> <p>Clock comes with remote control, used for setting up reminders.</p> <p>Reminders can be picture prompts & audio prompts.</p> <p>Large, clear screen</p> <p>Item is standalone (does not talk to other TEC items).</p> <p>Where you can find it on the Livity Life ordering system;</p> <p>Category: Reminders</p>

<p><u>Bellman - Vibrating Pad - B1270</u></p> <p><u>Item Code: CES65073-0053</u></p> 	<p>Device placed under users pillow and vibrates to wake upon linked alarm / alert (see other Bellman items in this catalogue section such as smoke alarms and push button transmitter).</p> <p>Vibrations are strong – should wake even heavy sleepers.</p> <p>Long 2m cable - easy to place where you want it</p> <p>This is not a standalone item – works in conjunction with other Bellman products. Does not talk to monitoring centre. Is not standalone – needs other Bellman products to talk to.</p> <p>Where you can find it on the Livity Life ordering system;</p> <p>Category: Living with Sight & Hearing Impairment</p>
<p><u>Push Button Transmitter B1420</u></p> <p><u>Item Code: CES65257-0053</u></p> 	<p>Item is water resistant. Can be used as a doorbell or person to person alert.</p> <p>Provides a vibration, flash to the person and can be used if the person is out of reach i.e., the garden where they wouldn't be able to hear their name being called</p> <p>This is not a standalone item – works in conjunction with other Bellman products, such as vibrating pad (CES65073-0053), pager (CES65308-0053) or visit flash (CES65356-0053).</p> <p>Does not talk to monitoring centre. Is not standalone – needs other Bellman products to talk to.</p> <p>Where you can find it on the Livity Life ordering system;</p> <p>Category: Living with Sight & Hearing Impairment</p>
<p><u>Bellman - Charger - B1260</u></p> <p><u>Item Code: CES65028-0053</u></p>	<p>Mains powered charger.</p> <p>Charges the pager (CES65308-0053) and pillow pad (Bellman Vibrating pad (CES65073-0053).</p> <p>This is not a standalone item – works in conjunction with other Bellman products,</p>

	<p>Does not talk to monitoring centre. Is not standalone – needs other Bellman products to talk to.</p> <p>Where you can find it on the Livity Life ordering system;</p> <p>Category: Living with Sight & Hearing Impairment</p>
<p><u>Hearing - Bellman Transmitter (B1430) B1431</u> <u>Item Code: CES65306-0053</u></p> 	<p>Transmitter connects to persons landline telephone and sends a signal to pager (CES65308-0053) when a call is received.</p> <p>This is not a standalone item – works in conjunction with other Bellman products, such as Bellman vibrating pad, pager (CES65308-0053) or visit flash (CES65356-0053).</p> <p>Does not talk to monitoring centre. Is not standalone – needs other Bellman products to talk to.</p> <p>Where you can find it on the Livity Life ordering system;</p> <p>Category: Living with Sight & Hearing Impairment</p>
<p><u>Hearing - Smoke Alarm - Ionisation - B1551</u> <u>Item Code: CES65314-0053</u></p> 	<p>Smoke detector transmitter for people who cannot hear an audio alarm. Audio alert works for other household members / visitors.</p> <p>This is not a standalone item – works in conjunction with other Bellman products, such as vibrating pad (CES65073-0053), pager (CES65308-0053), transmitter or visit flash (CES65356-0053).</p> <p>Alerts the person that the alert has been triggered via flashing lights and vibrations.</p> <p>Does not talk to monitoring centre. Is not standalone – needs other Bellman products to talk to.</p> <p>Where you can find it on the Livity Life ordering system;</p> <p>Category: Living with Sight & Hearing Impairment</p>

Bellman Visit Flash B1442**Item Code: CES65356-0053**

This Flash Receiver is a bright flashing light alerting the person that an alert has been triggered on a compatible Bellman item.

This is not a standalone item – works in conjunction with other Bellman products such as;

- Smoke alarm (CES65314-0053).
- Transmitter (often used as doorbell - (CES65257-0053).
- Telephone transmitter (CES65306-0053).

Vibrating pad or pager (CES65308-0053) so Person can still be alerted if asleep or out of room hosting the flash receiver.

Does not talk to monitoring centre. Is not standalone – needs other Bellman products to talk to.

Where you can find it on the Livity Life ordering system;

Category: Living with Sight & Hearing Impairment

Bellman Pager**CES65308-0053**

Portable pager alerting the person that an alert has been triggered on a compatible Bellman item.

This is not a standalone item – works in conjunction with other Bellman products such as:


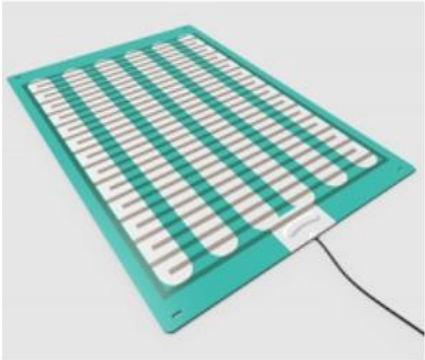
- Smoke alarm (CES65314-0053).
- Transmitter (often used as doorbell - (CES65257-0053).
- Telephone transmitter (CES65306-0053)



Vibrating pad or pager (CES65308-0053) so Person can still be alerted if asleep or out of room hosting the flash receiver.

This pager forms part of the Bellman paging system. The pager will vibrate and illuminate when activated. It is small enough to put in your pocket or clip to the side of your belt or skirt.


Different vibration patterns are produced for each alert, so blind and partially sighted people can also use this system. You can be notified to up to four different transmitters.

Does not talk to monitoring centre. Is not standalone – needs other Bellman products to talk to.

	<p>Where you can find it on the Livity Life ordering system;</p> <p>Category: Living with Sight & Hearing Impairment</p>
<p>General - Key Safe - (MS066) - C500 <u>CES15028-0098</u></p> 	<p>A small box, securely fitted to the outside of the property that holds the key.</p> <p>Enabling Carers / Responders to enter if a TEC alert is triggered.</p> <p>A key safe is fixed to an outside of the property using long screws.</p> <p>PRESCRIBER TO OBTAIN PERMISSION TO INSTALL BY HOME / PROPERTY OWNER.</p> <p>A key safe will only be agreed if ordered with TEC that requires monitoring or to support existing TEC in place.</p> <p>Where you can find it on the Livity Life ordering system;</p> <p>Category: Access</p>
<p><u>Enuresis sensor</u> <u>CES65749-1645</u></p> 	<p>Placed between bed and mattress, gives warning if moisture is detected.</p> <p>This is not a standalone item. Works in conjunction with Pager Buzz (approx. range 500 metres) (CES65748-1645).</p> <p>Carer informed via pager when moisture is detected.</p> <p>Can be wiped clean.</p> <p>Should be placed green and white side up.</p> <p>Does not talk to monitoring centre.</p> <p>Where you can find it on the Livity Life ordering system;</p> <p>Category: Detectors for urinary incontinence</p>
<p><u>Combined Occupancy/Enuresis) Sensor</u></p>	<p>Can be used as a combined occupancy (works on pressure) and enuresis sensor (moisture detector).</p>

<p><u>Item Code: CES65747-1645</u></p> 	<p>This is not a standalone item. Works in conjunction with Pager Buzz (approx. range 500 metres) (CES65748-1645).</p> <p>Carer informed via pager when moisture is detected, or pressure is no longer present in chair / bed.</p> <p>Where you can find it on the Livity Life ordering system;</p> <p>Category: Detectors for urinary incontinence</p>
<p><u>Move Sensor Voice Prompt CES 65663-0162</u></p> 	<p>A sensor which is triggered by movement.</p> <p>When sensor is triggered an audio message is played back.</p> <p>Family/ friend/ Carer can record one short message of approx. 20 seconds in length.</p> <p>Message will trigger and can act as a reminder / prompt a positive behaviour. For example, when next to the bed 'remember to use your walking stick', next to the front door 'Dad, don't forget your keys' etc.</p> <p>Standalone item.</p> <p>Cannot be linked to call monitoring centre.</p> <p>Where you can find it on the Livity Life ordering system;</p> <p>Category: Monitoring Movement in the home</p>
<p><u>Door & or window sensor – Linked</u></p> <p><u>Item Code: CES65730-0451</u></p>	<p>This is not a standalone item.</p> <p>The door/window alarm sensor can be set up in a number of different ways to alert / trigger when a door or window is opened or closed.</p>

	<p>Sensors can be set to trigger or not trigger at certain times. The prescriber needs to advise Liverty Life if they want timer's set-up and what the times are. For example, triggers of the door being opened may only be appropriate at night in some circumstances. The set up will depend upon the prescriber's assessment.</p> <p>The sensors will work with either;</p> <p>Lifeline Unit (CES65742-0451) – which will communicate to the Monitoring Centre. For monitoring purposes responders will be required who live within 30 mins of the property. (Monitoring charges may be applicable).</p> <p>Or</p> <p>Pager (CES65748-1645) – which will alert Carer if withing a 500-metre range.</p> <p>Where you can find it on the Livity Life ordering system;</p> <p>Category: Access</p>
<p><u>Telecare - Contact (Door / Window Sensor) - M132-105</u></p> <p><u>Item Code: CES65597-0162</u></p> 	<p>Kit comprises: 1 x pager, 1 x door/window transmitter with fixings.</p> <p>Alarm will trigger and a signal sent to the pager.</p> <p>Will alert carer if door or window is opened.</p> <p>Pager range is approx. 100 metres.</p> <p>Standalone item. Does not talk to the monitoring centre. If Pager (Carer) is outside of range the alert will not trigger.</p>

	<p>Where you can find it on the Livity Life ordering system;</p> <p>Category: Access</p>
<p><u>Floor Pressure Sensor Mat</u> <u>CES 65625-0280</u></p> 	<p>When pressure is applied to mat an alert will be triggered. Alert goes to pager.</p> <p>Mat is robust and hosts bevelled edges – to reduce trip hazard risk.</p> <p>Mat can be placed in various places as needed; next to bed, by the front door etc to alert Carer when the person may need support.</p> <p>The kit includes 1 x floor pressure mat sensor, 1 x alarm transmitter and a caregiver alarm pager.</p> <p>Standalone item (does not talk to any other TEC).</p> <p>Does not talk to the monitoring centre.</p> <p>Where you can find it on the Livity Life ordering system;</p> <p>Category: Monitoring Movement in the home</p>
<p><u>Oysta Pearl +</u></p> <p><u>Item Code: CES65762-1343</u></p>	<p>Device worn on person. Can be worn as a pendant alarm or in pocket.</p> <p>Comes with lanyard. Cannot be worn on the wrist.</p>



Device can be set to track GPS/ offers live location of where the person is.

Device can be set so alerts trigger if person travels outside specified area.

SOS button – Two-way voice communication

Falls detector – alerts responder if a fall is detected.

Reminders can be set– i.e., Medication reminders.

Prescriber to advise Livity Life on how they need this device set up. Prescribers please complete the Configuration form. This can be found:

https://www.centralbedsapp.co.uk/app-localpolicypractice/?hilite=matrix#Practice_Guidance

Prescriber of Pearl+ can contact IntegratingTEC@centralbedfordshire.gov.uk for support if required.

Where you can find it on the Livity Life ordering system;

Category: Monitoring Movement outside of the home